# APILUS 3G



The Gentle Way to 100% Permanent Hair Removal.

## APILUS 3G

### Clear Benefits

Even more effective, comfortable and versatile, the new Apilus 3G electrolysis devices will greatly impress you!

#### Gentleness and Security

Apilus 3G technology emits a perfectly stable 13.56 MHz radio frequency. It adapts to the variable impedance of the human body while delivering the appropriate amount of energy in order to achieve safe and gentle 100% permanent hair removal on all skin and hair types.

#### Performance, Precision and Customization

Apilus 3G treatment programs have been optimized for each body area and hair types, while taking into account the probe used. Your treatments will therefore be better adapted to your customers' needs. The calibration of the galvanic current and control of the units of lye have also been redesigned to provide an unparalleled efficiency during *Blend* treatments.

### **Exclusive Techniques**

#### MicroFlash

Rapid 13.56 MHz thermolysis technique working in milliseconds. Combined with insulated probes, it also allows the application of two proportional pulses (bulb, bulge). Comfortable, fast and effective, it is the technique of choice for all types of hair!




#### **MultiPlex**

Technique that combines slow thermolysis with a very short MicroFlash pulse. The slow thermolysis creates a light warming preparing the follicle to receive a rapid MicroFlash pulse. This intense action is ideal for the permanent destruction of difficult terminal hairs or even slightly distorted hairs.



#### **Omniblend**

Technique that combines galvanic current with a continuous high frequency current. The Apilus automatically adjusts the intensity of the currents depending on the sensitivity of the customer through the unique tolerance test function. It is comfortable and easy to use.

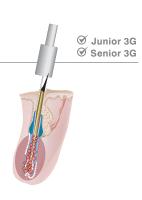

#### MultiBlend

Technique that combines the OmniBlend with one pulse of rapid thermolysis MicroFlash issued at the end of the process. An avantgarde modality specifically designed to permanently remove stubborn and distorted hair!



### 2 Versions to Suit your Needs

# APILUS | The Perfect | Introduction

The Apilus Junior 3G offers to the electrologist, in exemplary simplicity of use, the popular Flash and OmniBlend techniques as well as efficient and gentle treatments worthy of the Apilus brand. Undoubtedly, this introductory model will let you perform quality permanent hair removal in the most convenient way.

#### Modes

Thermolyse, Flash, Blend, OmniBlend, Electrolysis Anaphoresis, Cataphoresis



# APILUS | The Complete SENIOR 3G | Solution

The Apilus Senior 3G is the quickest, the most power ful and the most comfor table 13,56 MHz technology on the market. It integrates the exclusive iMM-Pac option and offers more than 10 treatment modalities as well as 1080 programs, allowing to adapt to all situations. Its high performance and versatility makes it the most popular electrolysis device in the world.

#### Modes

Flash, MicroFlash, MultiPlex, Blend, OmniBlend, MultiBlend, Electrolysis Anaphoresis, Cataphoresis



# APILUS 3G

### Used for over 30 years by thousands of centers around the world, Apilus is the number 1 brand in permanent hair removal.

	1	1
Features	JUNIOR	SENIOR
High performance circuit (25% more powerful than the Junior 3G)		•
Radio frequency 13.56 MHz	hundredths (sec)	thousandths (sec)
Superior comfort level		•
Treatment modalities : Thermolyse, Flash, Blend, OmniBlend, Electrolysis, Anaphoresis, Cataphoresis	•	•
Additional treatment modalities : MicroFlash, MultiPlex, MultiBlend		•
iMM-Pac technology including 2 exclusive modalities (MultiPlex and MultiBlend), a series of programs for insulated probes and a proportional pulsing function with variable delay.		•
Hair removal programs classified by modality, body area, hair size and probe type (insulated or non-insulated) allowing an unparalleled efficiency and ease of use.	512 programs	1080 programs
Possibility to create 60 customized programs		•
Cataphoresis programs for customized post-epilatory care	•	•
Unique tolerance test for more comfortable Blend treatments	•	•
Skin moisture sensor	•	•
"Processor Sentinel" protection circuit	•	•
UL key to calculate units of lye		•
Thermolysis "εL" energy level key		•
Automatic mode for operation without footswitch	•	•
Accessory test		•
Treatment management software with cumulative income statements		•
Insertion counter that compiles the number of treated hairs and cumulative treatment time	•	•
Probeholder connectors	1	2
Skin moisture sensor (not available in USA)	•	•

YOUR APILUS 3G DISTRIBUTOR

Made in Canada by: dectro INTERNATIONAL dectro.com

